



Cannulas & Trocars

Blunt Trocars


Used principally for the Hasson technique, the blunt trocar is packed with the Hasson fascia adaptor for primary port entry.

Features and benefits

- Smooth finish for easy attachment of fascia adaptor
- Blunt tip to prevent internal abdominal damage
- Lever for adjustable cannula depth
- Wide thread for port-site stability
- Used in patients with risk of adhesions from previous surgery
- Available in 10mm and 12mm diameters
- Tyvek packaging reduces space and waste per procedure
- Radiopaque markings ensure the distal tip shows up during X-ray

Blunt Trocars

Blunt trocar with smooth cannula and Hasson fascia adaptor for primary port entry.

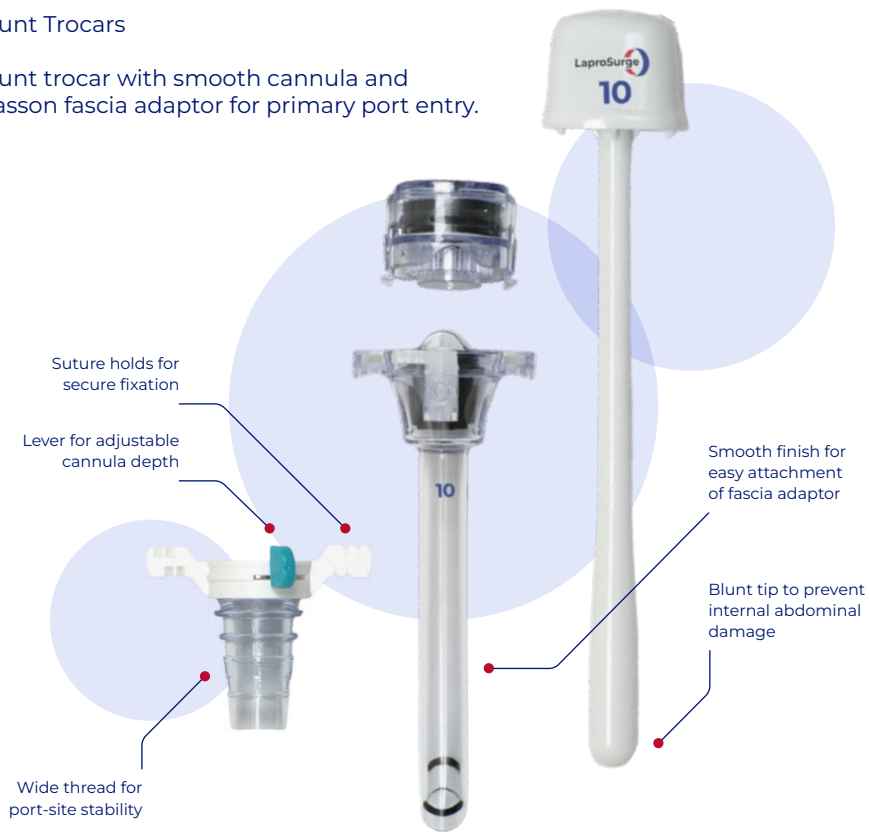

Product	Type	Code	Diameter	Cannula	Length	Pack	Box Quantity
Blunt Trocar	Standard	EC10 SHL AUTO	10mm	Smooth	100mm	Single	5
Blunt Trocar	Standard	EC12 SHL AUTO	12mm	Smooth	100mm	Single	5



+44 (0) 1923 220 711

info@laprosurge.com

LaproSurge Ltd, Unit 5b, Fishers Industrial Estate, Wiggenhall Road, Watford, Hertfordshire, WD18 0FN, England

UK Manufactured,
Globally Distributed,
Minimally Invasive.

laprosurge.com